

hadewe

Manual



SB-40

Directory

Foreword	3
Safety Instructions	4
Cleaning & Disinfection	8
Getting to Know the Unit	9
Before First Use	10
Area of Application	10
Operation	10
Application	12
Operation Steps	12
Working Positioning	13
Suction	15
Exchange Micro Filter	15
Exchange Control Unit Filter	16
Clamping the Bur/Tool	16
Setting up the Speed	16
Foot Switch	16
Maintenance & Care	17
Cleaning the Chuck	17
Maintenance Overview	18
Handling Instructions	19
Troubleshooting	20
Technical Details	21

Foreword

The SB-40 foot care unit was produced according to strict quality criteria and uses the latest technology in its design. The images in this manual may vary slightly from the product you have bought as the products are constantly being developed.

This user manual includes a detailed description and explanation of how to use the SB-40 unit. If you have any questions or suggestions, you can contact us by telephone or email.

This user manual belongs with the SB-40 unit. Keep it somewhere safe. If you give this product to someone else, please include this document as it includes important information on how to correctly use the product.

Short description of functions

The SB-40 is a device for the treatment of callus and nails. The handpiece is compatible with rotating tools such as diamond polishers, steel cutters and ceramic cutters. These can reach at speeds up to 40,000 rpm and therefore provide an efficient way of removing tissue. The direction of rotation can be changed by pressing a button. The integrated filtered suction function effectively sucks up any dust particles. The rotation speed and suction function can be adjusted using the keyboard. The device display will notify you when the filter needs to be changed and when the device is due for inspection.

Safety Instructions



The equipment should only be used by trained operators.



The installation of the working area has to correspond to the relevant regulations. Set up the unit in such a way, that the air slits are not blocked and the extracted air can escape well.



To avoid the risk of electrocution only connect the device to the mains power supply.



If you pull the power plug, it will stop the connection to the mains power supply. Set up the unit so that the power plug is always accessible. The power switch should also always be easily accessible.



Only switch on the unit if the handpiece rests secured in the handpiece holder, or if you hold it safely in your hand.



Prevent any liquids from entering into the control unit, since this could damage the electronic elements. Clean the outside of the box only with a damped cloth. Make sure that you have unplugged the unit before.



Only when pulling the mains cable, you can ensure that the device will be completely disconnected from mains. Check the mains cable regularly for damages and in case of such, it has to be exchanged or repaired by an authorized

Safety Instructions

service center.



Prevent your hair from wrapping itself around moving parts. If applicable, wear a hair net.



When working with materials, which might create dust or moisture, use a drilling unit with suction or spray technique. Wear a nose and mouth protection. Always follow the safety instructions, provided by the manufacturer of the materials being formed and used.



For safety reasons you must always wear disposable gloves and protective glasses when operating this unit. In addition, a face mask conforming to EN 14683 type IIR or EN 149 FFP 2 must also be worn.



Ensure that the patient's foot and the handpiece are held firmly in place during the operation. To prevent injury, work with extreme caution on the foot of the patient, and use only light pressure.



To ensure that the bur/tool sits firmly in the chuck, do not work with tools which:

- have an oily shaft,
- are worn,
- are bent.



Before using the equipment, assure that

- the tool can be pushed far enough into the

Safety Instructions

chuck (test this with a dipstick and clean if necessary),

- the inserted tool sits properly in the chuck and will not cant/jam while clamping it,
- the tool is approved for using it at the maximum speed (see manufacturer's instructions) or only use the speed required by the instrument,
- the tool is absolutely dry.



For safety reasons do not change the tool while the motor is running.



Switch off the unit immediately if

- the bur comes off during the treatment itself.
- the motor blocks during the treatment.



While doing any foot care, only operate the unit with the suction turned on.



Only use the suction function, if you have installed the suction filter. Check this filter every 10 operating hours. If it is full, please exchange it. A filter should be exchanged at least after 50 operating hours.



Never use the suction to vacuum any liquids (or air/water fog which contains particles of liquids).

Safety Instructions



Children should be supervised not to play with the appliance. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.

Always follow the operating instructions. Do not carry out any repairs, modifications or maintenance work yourself. This is only to be accomplished by an authorized specialist. The mains cable may also only be changed by the manufacturer, or an authorized specialist. Maintenance records may be supplied to repair workshops and trained technicians upon request. Never remove the housing, especially not while the equipment is plugged into the socket, since some components are under electric voltage. Any unauthorized or improper handling (e.g. trying to repair the unit) of the equipment will immediately nullify the guarantee/warranty and any other liabilities of the manufacturer.


Cleaning & Disinfection



The device is designed to be compatible with rotating instruments which do not penetrate the skin or mucosal or inadvertently may still cause injuries. For this application please adhere to the cleaning and disinfection instructions provided below. After every treatment please carry out the following steps:

- Remove the rotating tool (e.g. bur, capper etc.) first. Then handle the tool according to the instructions of its manufacturer.
- Remove tissue remains by wiping off the unit with a mildly moistened paper cloth. Afterwards clean all surfaces with disinfection cloth (e.g. Mikrozyd AF Tücher of Schülke and Mayr) There is no limit to the number of times this device can be cleaned and disinfected during its recommended lifespan.

Getting to Know the Unit

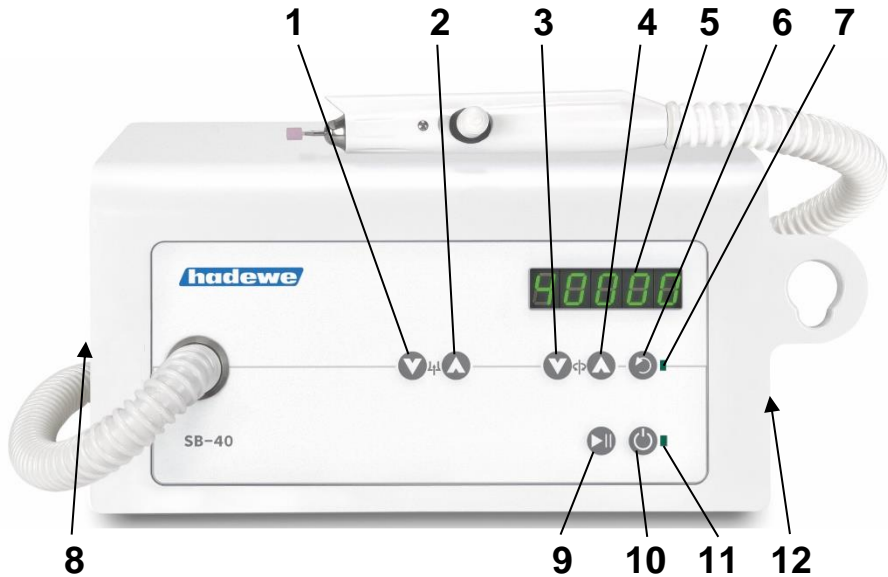

- 1 decrease the suction output
- 2 increase the suction output
- 3 reduce the speed
- 4 increase the speed
- 5 display of speed/suction level
- 6 change direction of rotation
- 7 LED is ON during reverse rotation
- 8 main switch
- 9 START/STOP switch
- 10 ON/OFF button (Standby)
- 11 standby LED
- 12 connecting socket for foot switch

Before First Use

Please read through the instructions carefully before using the unit and carefully note the safety advice. Always keep this manual accessible to any users of the drill.


Area of Application


This device should only be used by people who have completed a foot care qualification. Before use the operator should familiarize with the device by reading this manual or being instructed by another.

The appliance may only be operated in closed rooms. It can be used in a commercial or domestic environment.

The control unit can be set up in a cabinet or can be used as a portable unit (e.g. in a carry case). It is suitable for the following scopes:

- foot care (milling, trimming and polishing),
- manicure (grinding, milling and polishing works).

 The device is determined of operating rotating instruments, which will not penetrate the skin or mucous membrane if used appropriate, but could evoke damages accidentally.

 Caution: In the area of pedicure you have to work with suction techniques.

Operation

The working area must conform to the conditions described. Set up the unit in such a way, that the air slits are not blocked and the extracted air can escape well. To guarantee unhindered air flow there must a gap of at least 12 cm to the right of the device.

Position the device so that the operating buttons are easy to reach and the outlet of the hose is a maximum distance of 60 cm from the patient's foot. The unit should be 30 – 90 cm from the floor.

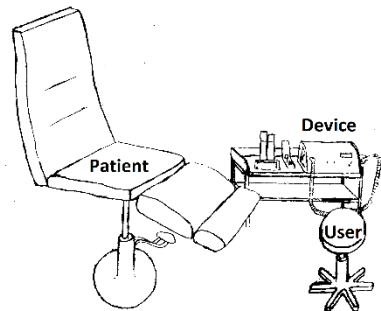


Image: Example of positioning of foot care unit in treatment room

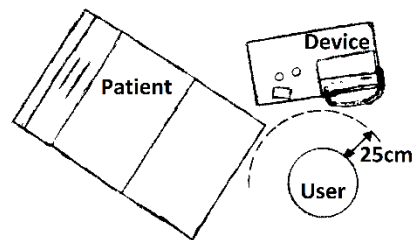



Image: Bird's eye view of work area

Before First Use

If you have a foot switch, insert its plug into the port on the right of the device which displays this symbol .

Check that the filter bag and dust filter are correctly fitted. Ensure that the tool is fixed.

On the left side in the back of the device you will find the main power switch. Press this switch to turn the unit on.



The device is now on standby. Pick up the handpiece. Should the handpiece be twisted by the hose you can adjust that at the control unit.



To use the handpiece and suction function press the button (9) on the front of the unit.

Application

Operation Steps

In the following section the order in which you should use the device is described. Additionally, the individual steps will be explained in detail. Read the complete manual before use.

After having clamped the instrument, press the standby key (10) to switch on the device. Use the buttons (3) and (4) to select a rotation speed. Manufacturers often provide recommendations for optimal rotation speeds. However, it is important that you never exceed the maximum recommended rotation speed. Now press the start/stop-button (9) to set the instrument in rotation.

You can now choose a rotation direction using the button (6). The most commonly used direction is clockwise. Most milling instrument (steel instruments with steel cutters) are optimised for clockwise use. The rotation direction is irrelevant for diamond cutters and high-grade corundum sanders in terms of efficiency. Working anti-clockwise can be beneficial if, for example, the dust obscures your field of vision when working clockwise.

Using buttons (1) and (2) you can adjust the intensity of the suction. After choosing a rotation speed and suction power you can start working.

If you want to take a break press the start/stop-button (9). The handpiece motor and suction will turn off. The handpiece can now be placed in its

holder on the right of the unit. If you want to end the work completely, use the standby key (10). The display will then also switch off.


Alternatively, you can switch the device to standby by pressing the button on the foot switch.

We recommend you disconnect the device from the mains power supply if it is not going to be used for a while. This will save electricity. The mains plug must also be removed when it comes to retiring the device.

Application

Work Positioning

The following images are examples of how to work on various parts of the foot. Always make sure that the foot is firmly gripped and that you can work freely and without strain.

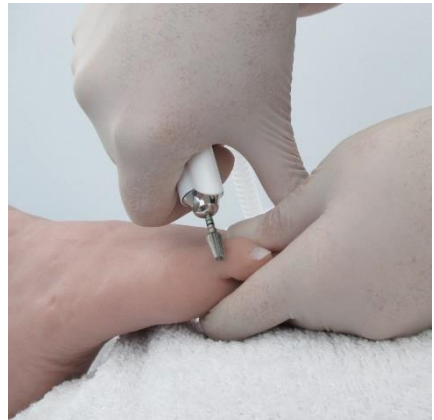

Application



Application

Suction

Switch on the suction with button (2). The suction output can be regulated in 5 levels. You will see the selected level in the display. Level 0 would mean that the suction is turned off.

In case the suction motor overheats, the unit will switch off automatically. Leave the unit to cool off for a while before starting it again. Should the same failure appear again, please return the unit to your supplier for inspection.

Exchange Micro Filter




When changing the filter, it is advisable to wear a face mask conforming to EN 14683 type IIR or EN 149 FFP 2.

Only operate the unit with a suction filter assembled. Control the filter bag after every 10 operating hours. Replace a full filter bag immediately, but no later than **FC** appears on the display. You will receive this message every 50 hours of operation in the suction mode.

You can read the operating time since the last filter change by pressing and holding button (6) while pressing the power button at the same time. First, the five-digit number of the software version of your device appears, then the operating time since the last filter change is displayed.

At first, unfold the paper of the paper

filter (art. 5115) bag and push your finger into the opening in the front.

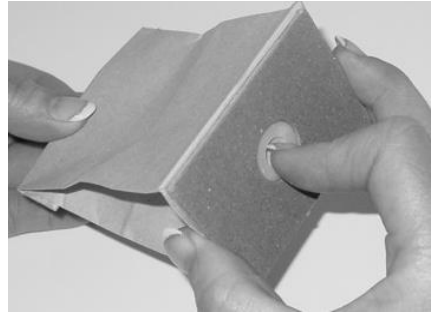

Pull out the lid of the dust chamber and remove the old filter.



Hold the filter in the middle to prevent it from bending and sliding off.




As an option you may also use the micro filter (art. 5179) instead of paper filter (art. 5115).

Application

After the filter bag has been changed, the operating time must be set to zero. To do this, set the device to standby by pressing button (10). Now hold down button (1) and (2) until **FC** appears in the display.


Exchange Control Unit Filter

This filter (art. 3752) has to be exchanged after every 200 operating hours. Remove the old filter and place the new one in its holder (see picture).

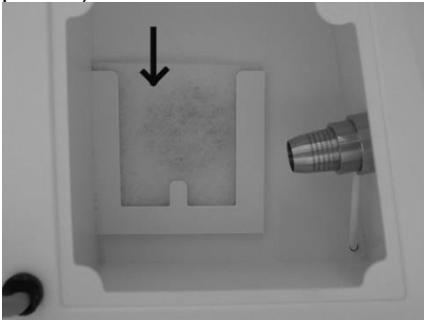

Clamping the Bur/Tool

The handpiece is equipped with a chuck, which can hold all burs with a diameter of 2.35 mm (this is the diameter usually used for tools in podiatry). To open the chuck, push the knob forward with your thumb and leave it there. Now insert the tool as far as possible. Pull back the knob and the tool will be firmly clamped.

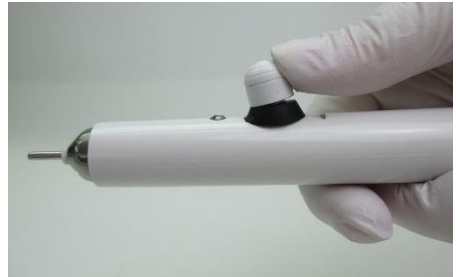


For security reasons, only change the tool, if the handpiece is turned off.

Setting up the Speed

Hold the handpiece (incl. tool) in your hand or store it in the handpiece holder attached to the side of the control unit. Switch on the unit with button (9). Now you can also set the speed with buttons (1) and (2) and choose its direction of rotation with button (6).

Foot Switch (art. 0970)

Operation of the foot switch

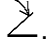Only use the original hadewe foot switch, to avoid any damages. Insert the plug of the foot switch into the socket with the symbol .

If you press the foot switch for a longer period, you will switch on/off the complete control unit.

Maintenance & Care

Cleaning the Chuck

In the process of time, some dirt might collect itself within the hand-piece, and might affect its functions. Therefore, the chuck has to be cleaned monthly.

How to remove the chuck

Screw off the top of the hand-piece with the spanner shown in the picture:




Now clean the inside of the hand-piece and the top with the brush.

For opening the chuck, push the button forward:



Set the screw spanner (wrench) on the axle and hold it tight. Now use the front of the second spanner to unscrew the chuck:

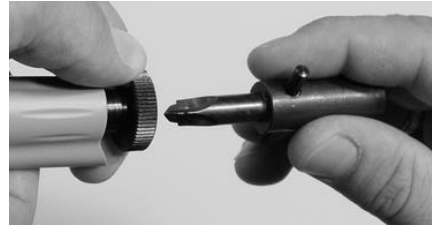

To assemble the chuck, follow the steps in reverse order. Make sure that the chuck is opened during the assembly (the handpiece button must be in forward position) and is firmly applied afterwards.

How to clean the chuck

Soak the chuck in ethyl alcohol or a cleaning liquid for tools. Clean the slits of the chuck with the brush:



Use the brush to remove all dirt residues from the drill-hole of the chuck:



Rinse the chuck very well, and make

Maintenance & Care

sure you let it dry well afterwards. Cover your index finger with a very little amount (just a drop) of oil and rub in between your thumb and index finger. Dab off the oil until your fingers is only covered with a very thin film of oil. Now move the chuck in between your thumb and index finger. There should be only a very thin film of oil left on the chuck. Please act according to the saying: Less is more.

Maintenance Overview

What has to be done?	How often? Done by whom?
cleaning the chuck	monthly by user
changing the filter 5115 or 5179	every 50 hours of operation by user (check after each 10 hours)
changing the control unit filter	every 200 hours of operation by user
maintenance by re- pair workshop (changing worn parts; cleaning; check leakage cur- rent)	every 250 hours of operation, af- ter one year at the latest by re- pair workshop

Handling Instructions



Never put the handpiece into any kind of liquid, and never oil or grease the handpiece. Clean its surface with a dampened cloth. Pay attention, that no moisture or liquid enters the handpiece. Never use spray disinfection.



Never put any wet burs or instruments into the handpiece. Its moisture could damage the interior parts of the handpiece.



Any materials used by hadewe prevent a formation of rust. In cases in which a rust infestation was found within the handpiece, thorough inspections and examinations always have proved that the rust had a foreign origin (which was usually caused by the tools used in these circumstances). Even in high-grade steel instruments of well known manufacturers, formation of rust can appear, e.g., through wrong treatment with chemical cleaning and disinfection liquids. This extraneous rust affects internal parts of the handpiece in form of rust bloom and causes pitting. Therefore, when you not working with the handpiece, always remove the bur/tool from it.



Always make sure that you do not let the handpiece fall down or exposes it to any shocks. The ballbearings and the motor could be damaged by this. Only work with tools, which are flawless. Do not work with tools

which

- are bent or out of balance
- whose shaft is worn out or
- are rusty.



Working with faulty tools will cause a stronger vibration of the handpiece. A damage of the ball bearings, clamping technique or motor could be the consequence.



Working with high pressure will not improve the outcome of your work, since the motor will even be slowed down. A second disadvantage is the stronger impact on the ball bearings and the motor, which will shorten their life span in the long run.



In case of sending the device by mail remove the dust filter before.



Only suck up dust or milling waste, never suck up larger particles such as small bits of cotton wool as this may be get caught in the handpiece and lead to blockages.

Troubleshooting

Handpiece vibrates, is very noisy and/or becomes warm in the front part.

Use of bent burs/tools. → Use a different bur/ tool to test this. The bur/tool itself can be tested with the hadewe Bur Shaft Tester (art. 4990) to see if it is bent.

Maximum allowed speed of the tool was exceeded, which caused a stronger vibration. → Pay attention to the tool (e.g. bur, capper and instrument) manufacturers handling details.

Bur/tool can be inserted only with difficulty or not at all.

Use of bent burs/tools. → Use a different bur/tool to test this. The bur/tool itself can be tested with the hadewe Bur Shaft Tester (art. 4990) to see if it is bent.

Chuck is dirty. → Clean the chuck according to the instructions in the manual.

Bur does not clamp or does not spin.

Chuck is not tightened well. → Tighten the chuck firmly according to the instruction in the manual.

Knob of the handpiece cannot be moved easily.

Chuck is dirty. → Clean the chuck according to the instructions in the manual.

Unit does not work at all.

Check, if the main switch of the unit is turned on at the side of the unit.
Foot switch is used constantly or faulty.
→ Pull out the socket of the foot switch to see if this is the cause.

Pull out the mains plug to reset the device.

The unit shows different sound levels at different sites of use.

The sound level varies according to the surroundings of the unit. If you set up the unit directly in front of a wall, the sound is reflected stronger, as if it would, while standing further away from a wall.

Handpiece motor stutters shortly (approx. 1 sec) when turning unit on.

In very rare cases this effect did appear. It will not harm or damage the device.

Suction function doesn't work.

The suction motor will turn off automatically if it overheats. Let the device cool down before using it again. Send the device off for inspection if this occurs repeatedly.

Message "FC" on the display

Changing the filter (Filter Change) → Page 16

Error "E1" on the display

The speed of the handpiece motor deviates too much from the set value.

Error "E2" on the display

The speed of the suction motor deviates too much from the set value.

Error "E3" on the display

The motor current from the handpiece has been too high for too long.

Technical Details

0987 SB-40

total weight: 3,8 kg
dimension: W273xH142xD186 mm
incl. exterior parts: W313xH142xD235mm
input: 230 V~, 50 Hz
power consumption: 0,9 A
vacuum: 2,9 - 5,8 kPa
fuse primary:
 F1: 315 mA time delayed safety fuse
 F2, F3: 3,15 A time delayed safety fuse
fuse sec.: thermal fuse
protection class II
This device contains a ground connection for functional purposes only
accuracy of speed display: $\pm 10\%$

handpiece:
type B applied part
diameter: 17-22 mm
length: 142 mm
speed: 6.000 - 40.000 rpm
vibration $< 2,5 \text{m/s}^2$

Accessories & Spare Parts

art.	
5115	paper filter
5179	micro filter
3752	filter
0970	foot switch
6158	cleaning brush

The clamping system is designed for rotating instruments, which have a shaft of 2,35 mm in accordance with DIN EN ISO 1797-1

Technical Details

Ambient Conditions

The device is intended for use in dry closed rooms.

Operation:

Temp.: +10°C to +26°C
Rel. humidity: 25% to 75%
(not condensed)
Air pressure: 700 hPa to 1060 hPa

Storage:

Temp: -5°C to +55°C
Rel. humidity: 10% to 95%
(not condensed)
Air pressure: 500 hPa to 1060 hPa

Transport (up to 4 weeks):

Temp: -5°C to +55°C
Rel. humidity: 10% to 95%
(not condensed)
Air pressure: 500 hPa to 1060 hPa

Operating Time

The device is designed for extended usage.

Declaration of Conformity

We declare under our sole responsibility, that the product (0987) podiatry unit, to which this declaration relates, is conform to the following standards or normative documents:

2014/30/EU directive EMC
2014/35/EU low voltage directive

hadewe GmbH – Grambartstraße 10 –
30165 Hannover

Accessories & Spare Parts

art.

5179	micro filter
3752	control unit filter
0970	foot switch
6158	cleaning brush
5179	micro filter


The clamping system is designed for rotating instruments, which have a shaft of 2.35 mm in accordance with DIN EN ISO 1797-1.

Safety Check

A safety check must be done in compliance with the regulations of the country where the device is used. The inspection interval is once a year.

Technical Details

Disposal

 The device and filter may contain infectious material such as residual sanding dust or skin particles. Therefore, please bear in mind the following disposal instructions.

Disposing of the device

Old devices bought in Germany that correspond to the electronic device law should be sent directly to the manufacturer (hadewe). We will dispose of the device free of charge. The devices should not be taken to public disposal plants (WEEE Reg.-Nr. DE20392713, b2b device). For disposal of this unit out of Germany, please contact the site where you bought this hadewe product.

Note for the disposal of the device/ treatment system

Wear protecting gloves of the class 2 in accordance with DIN 374-2 during the disassembling of the device. Also, wear a protective mask for your face according EN 14683 version IIR or EN 149 FFP2.

Disposal of dust filter and coarse dust filter

The filter bag and the coarse dust filter are considered low contaminated waste. Dispose those in in transparent, liquid- tight, resistant plastic sacks and do not dispose them compressed in the regular household garbage.

Symbols



Caution!



Important Note!



Follow instructions for use!



temperature limit
(indicates the lower and upper temperature limit)



air pressure: permissible range



air humidity: permissible range



Keep dry!



ON/OFF button (standby) (no disconnection from main power supply)



START/STOP button



change direction



suction



rotation



alternating current



direct current



safety fuse



socket foot switch



class II



date of manufacture, name and address of the manufacturer



Do not dispose in household waste

hadewe

hadewe GmbH
Grambartstraße 10
30165 Hannover
Germany
info@hadewe.de
www.hadewe.de